

# Organizing your image files

---

## Overall folder organization


---

- Parent folder
  - Experimental\_condition\_1
    - Image\_replicate\_001\_channel1.TIF
    - Image\_replicate\_001\_channel2.TIF
    - Image\_replicate\_002\_channel1.TIF
    - Image\_replicate\_002\_channel2.TIF
    - ...
    - Image\_replicate\_015\_channel2.TIF
  - Experimental\_condition\_2
  - ...
  - Experimental\_condition\_N
  - data\_information.xlsx

## Single folder = 1 unique experimental condition

---

- A single folder should contain all replicates for a given experimental condition.
- Replicates should have the following file syntax:
  - (meta-data)\_(unique replicate)\_(w)\_(wave-length information)
  - Example of 3 replicates imaged across 3 channels with unique IDs **002** to **004**.

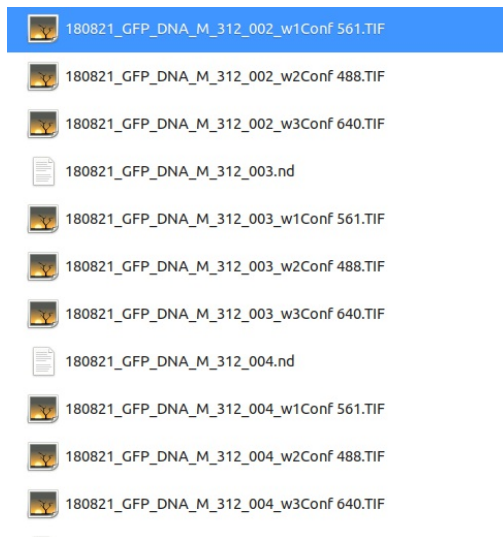

## Single parameter sweep = List of folders

---

All the different experimental conditions should be under a parent folder with an **excel file** outlying important data.

- Example organization of a MED1-IDR concentration sweep shown below:



- data\_information.xlsx contains 4 columns:
  - folder\_path - Lists all relevant folders (see previous image)
  - Channel\_number - Gives name to channel, not required to be unique. If a file doesn't exist in a channel, leave the cell blank.

### **Examples**

	A	B	C	D	
1	folder_path	Channel 488	Channel 561	Channel 640	
2	GFP_DNA20_MED1_10	GFP_M10	MED1_10	DNA20	
3	GFP_DNA20_MED1_19	GFP_M19	MED1_19	DNA20	
4	GFP_DNA20_MED1_39	GFP_M39	MED1_39	DNA20	
5	GFP_DNA20_MED1_78	GFP_M78	MED1_78	DNA20	
6	GFP_DNA20_MED1_156	GFP_M156	MED1_156	DNA20	
7	GFP_DNA20_MED1_312	GFP_M312	MED1_312	DNA20	
8	GFP_DNA20_MED1_625	GFP_M625	MED1_625	DNA20	
9	GFP_DNA20_MED1_1250	GFP_M1250	MED1_1250	DNA20	
0					
1					